

Zeta Plus™ C(P) & S(P) Series

Lenticular Depth Filter Cartridges

Zeta Plus C, CP, S & SP Series filters are a family of advanced depth filters designed to provide optimal clarification of pharmaceutical, biological, bioprocess, cosmetic and food & beverage fluids. Zeta Plus C, CP, S & SP Series filters are composed of inorganic filter aid, cellulose and a resin system that imparts a positive charge to the filter matrix, aiding its adsorption of negatively charged contaminants.

Grade Choice


The Zeta Plus C, CP, S & SP Series filter are available in several grades, that meet the different filtration needs required by your applications. Filtration ratings are given in the following table for the various grades. The optimal filtration system for your particular application can be determined by on-site testing or by sample evaluation by a member of our Scientific Application Support Services group (SASS).

Pharmaceutical CP and SP Grade

Zeta Plus CP and SP media are manufactured following procedures described in 3M Purification's Drug Master Files submitted to the FDA (United States Food and Drug Administration). The characteristics and filtration mechanisms of the CP and SP Grades are the same as the standard C and S Series. The CP and SP Grade has complete product control and traceability.



Picture 1 : Zeta Plus Family

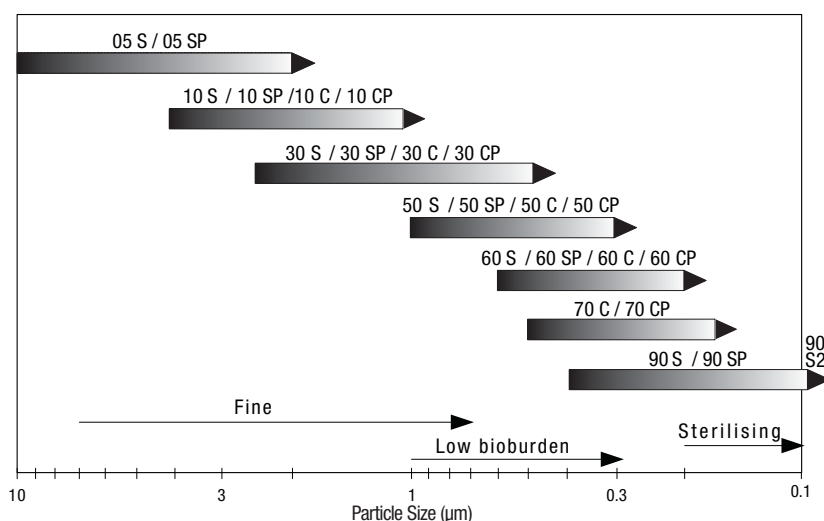

Graphic 1: Zeta Plus C, CP, S and SP Grade Selection

Applications

Zeta Plus C, CP, S & SP Series filters are optimised for clarification of pharmaceutical, biological, bioprocess, cosmetic and food & beverage fluids. Examples of applications include:

Bio Pharmaceuticals

- Blood Products
- Antibiotics
- Cough Syrups
- Intravenous Solutions
- Growth Media
- Bioprocess Fluid
- Cosmetics

Food & Beverages

- Spirits
- Water

Features and Benefits

Depth filtration media designed to retain contaminants by mechanical entrapment and electrokinetic adsorption

- High contaminant holding capacity for economical filtration and reliable particle removal

Full range of scaleable capsule and cartridge filter configurations

- Allows pilot testing and scale-up with the same materials that will be used in full-scale systems

FDA Drug Master File and USP Class VI Biological Safety

- Eases validation and regulatory submissions by providing vital documentation and traceability

Self-contained, single-use disposable capsule and cartridge modules

- Reduced labour time for change-outs and capsules eliminate cleaning validation



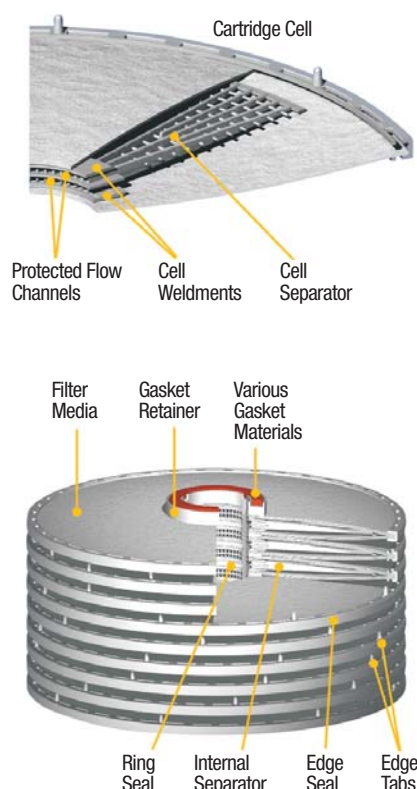


Figure 1: Zeta Plus Cartridge Construction



During the manufacturing process, molecules carrying a positive charge are chemically bonded to the filter matrix permanently forming an interconnected filtration structure with positively charged electrokinetic capture sites. The resulting porous depth filter is a tortuous network of charge-enhanced flow channels capable of removing bacteria, particulate, cellular debris and submicronic contaminants to a level which mechanical screening alone cannot achieve.

Cartridge Construction


Zeta Plus C, CP, S and SP Series cartridges are designed for use with 3M Purification Zeta Plus sanitary style stainless steel filter housings. Cartridges are constructed from individual cells of Zeta Plus C, CP, S and SP Series filter medium assembled together with polypropylene ring seals under predetermined compression and unitised by three 316 stainless steel bands. Each cell is constructed using polypropylene molded edge seals and separators for high performance. Various gasket materials are available depending upon application. Filter cartridges are available in 8, 12 and 16-inch diameters.

Capsule Construction

Zeta Plus C, CP, S and SP Series filters are available in ready to use, fully disposable Zeta Plus BC capsules. The Zeta Plus BC Series capsules are completely encapsulated disposable filters for small volume biological, bioprocess and pharmaceutical filtration. Zeta Plus BC capsules are autoclavable and available in three sizes (25, 650 and 1300 cm² filtration area) to provide a high degree of application flexibility. They provide a perfect tool for scaling experiments.

Efficient Flow Rates

The graph shows initial Δp values for Zeta Plus CP and SP Series filter media versus specific flow rates (litre per min/m²), obtained with clean water.

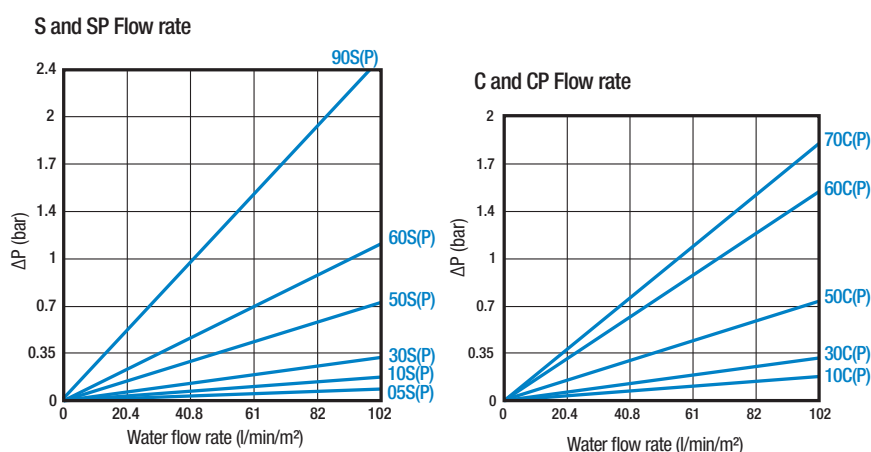

Graphic 2: Zeta Plus Flow Rate Graphs

Quality Control

The Zeta Plus C, CP, S & SP Series cartridges and sheets are manufactured following an ISO quality certified assurance scheme. All the cartridges are labelled with a lot identification number to provide complete traceability.

Sanitary Housings

3M Purification provides a wide range of standard and specially designed sanitary filter housings to hold Zeta Plus C, CP, S and SP Series filter cartridges. All housings are designed with the pharmaceutical and bioprocessing industries in mind and feature 316L mirror-polished surfaces and easy-to-clean components.

Regulatory Support

Zeta Plus C, CP, S & SP Series filters are produced from materials which are 21CFR listed and filter components have successfully passed USP Class VI Biological Reactivity Tests.

Compatibility

Temperature, contact time, chemical concentration, flow rate, aeration and differential system pressure are all variables that can affect chemical compatibility of filter materials. Therefore it is always advisable to perform tests replicating process conditions to determine ultimate compatibility. Table 1 is provided as a guide only. Please contact your local 3M Purification SASS Team for further support on validation and compatibility.

Table 1: Chemical Compatibility						
CHEMICALS	S(P) Max. (82°C)	C(P) Max. (60°C)	GASKETS			
			Nitrile	Fluorocarbene	Ethylene Propylene	Silicone
Ethyl Acetate	G	G	P	P	P	F
Acetone 100%	G	G	P	G	G	G
Acetic Acid 5%-20%	G	G	G	G	G	G
Benzene	G	G	P	F	P	F
Chloroform	G	G	P	G	P	F
Chlorothene	G	G	P	G	P	F
Methylene chloride	G	G	P	F	P	P
1,4 - Dioxane	G	G	P	P	P	P
DMF	G	G	F	P	P	No Data
Water (25°C - 82°C)	G	G	G	G	G	G
Ethanol 10%	G	G	G	F	G	G
Ethanol 50%	G	G	G	F	G	G
Hexane	G	G	G	G	P	F
Sodium Hydroxide 2%	F-P	F-P	F	G	G	G
Methanol	G	G	G	P	G	G
Methyl Ethyl Ketone (MEK)	G	G	P	P	G	G
Methyl Isobutyl Ketone (MIBK)	G	G	P	P	G	G
n-Heptane	G	G	G	G	P	F
Carbon Tetrachloride pure	G	G	F	G	P	F
1,1,1 - Trichloroethane	G	G	P	G	P	F
Xylene	G	G	P	G	P	F

Explanation of ratings:
 G = Satisfactory - to maximum temperature given for material, unless restriction is noted.
 F = Fair - to maximum temperature given for material, unless restriction is noted.
 P = Not recommended.

Extractables

A static soak test of Zeta Plus C, CP, S & SP materials was performed to determine the maximum amount of organic and inorganic compounds extracted by various solvents. Extraction tests were performed using 100 ml of fluid per 20 cm² of media.

Operating Parameters

Table 3: Recommended Operating Parameters					
	Maximum Operating Pressure	Maximum Operating Temperature	Recommended Flush	Sterilization Parameters	Maximum Flow Rate (FLUX) (in l/min per m ²)
Standard and Special precoat Zeta Plus Cartridges	2.4 bar	82 °C	54 l/m ² at 20 litre per min/m ²	autoclave or in-situ steam sterilisation for 1 hour at 121 °C	1.2 to 12
Zeta Plus BC25 Capsules	Inlet pressure: 2.75 bar Media pressure: 2.4 bar	40 °C	54 l/m ² at 20 litre per min/m ²	Autoclave for 30 minutes at 121 °C (1 cycle)	1.2 to 12
Zeta Plus BC1000 and BC2000 Series Capsules	Inlet pressure: 5.5 bar at 25 °C 2.75 bar at 60 °C Media pressure: 2.4 bar	60 °C	54 l/m ² at 20 litre per min/m ²	Autoclave for 30 minutes at 121 °C (up to 3 cycles)	1.2 to 12

Table 2: Zeta Plus C(P) & S(P) Series extractables				
Zeta Plus S(P) Series extractables (mg/g of media)				
Solvents	Na	K	Ca	Mg
Water	0.075	0.039	0.160	0.070
Water (82°C)	0.077	0.133	0.066	0.014
Ethanol 10%	0.013	0.004	0.017	0.002
Ethanol 50%	0.014	0.002	0.012	0.002
Acetic Acid	0.062	0.157	0.490	0.936
n-Heptane	<0.001	<0.001	<0.001	<0.001
Solvents	Fe	Heavy metals		Organics
Water	0.0009	<0.0001		<0.005
Water (82°C)	0.0011	<0.0001		<0.005
Ethanol 10%	<0.0005	<0.0001		<0.005
Ethanol 50%	<0.0005	<0.0001		<0.005
Acetic Acid	0.170	<0.0001		<0.005
n-Heptane	<0.0005	<0.0001		<0.005
Zeta Plus C(P) Series extractables (mg/g of media)				
Solvents	Na	K	Ca	Mg
Water	0.175	0.042	0.045	0.008
Ethanol 10%	0.210	0.065	0.080	0.014
Ethanol 50%	0.165	0.025	0.063	0.010
Acetic Acid	0.255	0.130	0.115	0.019
n-Heptane	0.007	0.003	0.012	0.002
Solvents	Fe	Heavy metals		Organics
Water	<0.0005	<0.0001		0.031
Ethanol 10%	0.0011	<0.0001		0.014
Ethanol 50%	0.0008	<0.0001		0.009
Acetic Acid	0.0050	<0.0001		0.067
n-Heptane	<0.0005	<0.0001		0.002
< = limits of detection				
NOTE : Extractable levels will vary from sample to sample as the raw materials are natural products.				

Zeta Plus C, CP, S and SP Series Ordering Guide

1. Cartridges

Zeta Plus	Diameter (inch)	Cartridge Design	Gasket Type	Grade	Formulation
Z	08	P = Plug-in 7 cells (0.23 m²) P2 = Plug-in 2 cells (0.065 m²) P4 = Plug-in 4 cells (0.13 m²)	Standard A = Silicone (MVQ)*	05	S SP
		D = Standard 8 cells (0.26 m²)	Standard D = Nitrile (NBR)*	10 30 50 60	S SP C CP
	12	C = 9 cells / small (0.85 m²) B = Special precoat 12 cells (1.1 m²) D = Standard 16 cells (1.5 m²) S = Special 7 cells (0.7 m²)	Options A = Silicone (MVQ)* B = Fluorocarbon (FPM)* C = EPR (EPDM)*	70	C CP
	16	M = Standard diffusion netting 14 cells (3.2 m²) D = Standard 15 cells (3.4 m²) S = Special precoat 9 cells (2.1 m²)	*ISO Designation	90	S SP

2. Zeta Plus BC Capsules

Catalogue Number	Grade and Formulation																								
BC 0025L (Luer) - 25 cm² area BC 0025S (Sanitary) - 25 cm² area BC 1000A (single pack) - 650 cm² area BC 2000A (single pack) - 1300 cm² area	<table><tr><td>10 C</td><td>10 CP</td><td>05 S</td><td>05 SP</td></tr><tr><td>30 C</td><td>30 CP</td><td>10 S</td><td>10 SP</td></tr><tr><td>50 C</td><td>50 CP</td><td>30 S</td><td>30 SP</td></tr><tr><td>60 C</td><td>60 CP</td><td>50 S</td><td>50 SP</td></tr><tr><td>70 C</td><td>70 CP</td><td>60 S</td><td>60 SP</td></tr><tr><td></td><td></td><td>90 S</td><td>90 SP</td></tr></table>	10 C	10 CP	05 S	05 SP	30 C	30 CP	10 S	10 SP	50 C	50 CP	30 S	30 SP	60 C	60 CP	50 S	50 SP	70 C	70 CP	60 S	60 SP			90 S	90 SP
10 C	10 CP	05 S	05 SP																						
30 C	30 CP	10 S	10 SP																						
50 C	50 CP	30 S	30 SP																						
60 C	60 CP	50 S	50 SP																						
70 C	70 CP	60 S	60 SP																						
		90 S	90 SP																						

3. Sheets

By request, Zeta Plus media can also be supplied in sheet format. Please contact your local 3M Purification representative for further details on how to order.

Important Notice

3M Purification MAKES NO WARRANTIES, EXPRESS OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. Since a variety of factors can affect the use and performance of a 3M Purification product in a particular application, some of which are uniquely within the user's knowledge and control, user is responsible for determining whether or not the 3M Purification product is fit for a particular purpose and suitable for user's method of application.

Limited Warranty

3M Purification warrants it this product to be free from defects in material and workmanship during normal use for a period of one (1) year from the date of shipment from the factory. If the Product(s) is (are) defective within this warranty period, your exclusive remedy and 3M Purification's sole obligations shall be, at 3M Purification's option, to replace or repair the Product(s) or refund the original purchase price of the Product(s). This warranty does not apply to failures that result from abuse, misuse, alternation or damage not caused by 3M Purification or failure to properly follow installation and use instructions.

Limitation of Liability: 3M Purification will not be liable for any loss or damage arising from the use of the Product(s), whether direct, indirect, special, incidental, or consequential, regardless of the legal theory asserted, including warranty, contract, negligence or strict liability. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights and you may have other rights which vary from state to state, or country to country.

Zeta Plus and 3M are trademarks of 3M used under license.



3M United Kingdom plc

3M Purification

3M Centre, Cain Road
Bracknell, RG12 8HT
United Kingdom
Tel.: 0845 602 5237
Fax: +44 (0)1344 858559
E-mail: innovation.uk@mmm.com

Your specialist dealer:



FUHR GmbH
Am Weinkastell 14
D - 55270 Klein-Winternheim
Tel. 06136/9943-0
Fax 06136/9943-25
E-mail: info@fuhr-gmbh.com
www.fuhr-gmbh.com/3m.htm

For more contact addresses visit our website www.3M.eu/filtration or www.3Mpurification.com/international.

Data may be subject to change without further notice.

© 3M 2010. All rights reserved.

v.2 - DOC10021
LITCTZPCS1.EU - 1110